

Antelope Valley Union High School District
TwinJect Autoinjector Epinephrine Injection Device



FIRST DOSE
ADMINISTRATION

Remove caps labeled
"1" and "2."


Place rounded tip against the
middle of the outer thigh. The
needle may be injected through
clothing. Press down hard until
auto-injector fires. Hold in place
while slowly counting to 10.

Immediately **call 911**.

Never place thumbs, fingers,
or hands over rounded tip.

SECOND DOSE
ADMINISTRATION

If symptoms have not improved
within about 10 minutes after
the first injection, a second
dose is needed.

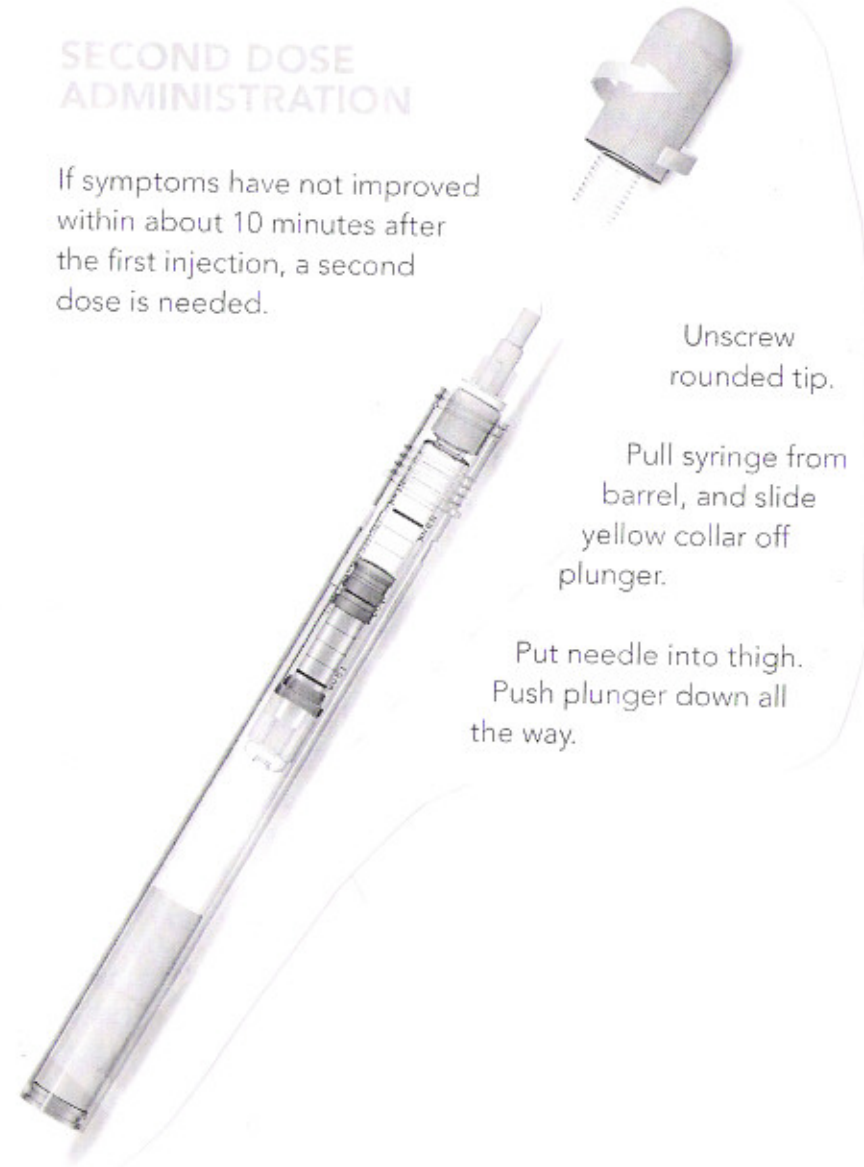

Unscrew
rounded tip.

Pull syringe from
barrel, and slide
yellow collar off
plunger.

Put needle into thigh.
Push plunger down all
the way.